



Cure Medical supports medical research. For information on Cure Medical catheters, visit www.curemedical.com or contact your distributor of quality healthcare products.



Instruction Guide to Clean Self-Intermittent Catheterization For Men Using the Cure Catheter® with Coudé Tip


This educational material is provided by Cure Medical in an effort to answer questions and address possible concerns about the intermittent catheterization process. By following the simple steps for clean self-intermittent catheterization, it is possible to quickly master the technique and make the process part of your daily routine for bladder management and urinary tract health.

Cure Medical wishes to thank Anne Boisclair-Fahey, DNP, RN, CPNP for writing this educational material. Her experience as a Senior Director of Advanced Practice Providers (APPS) at Children's Minnesota hospital, combined with over 20 years of clinical expertise in pediatric urology as a pediatric nurse practitioner, make her contribution to this instruction guide invaluable.

This instruction guide is not a substitute for medical advice from your healthcare provider.

The Urinary System

The urinary system contains two kidneys, two ureters, the bladder and the urethra. The kidneys filter the blood and produce urine. The urine travels from the kidneys down the ureters and into the bladder, where it is stored until emptied during urination or catheterization. The urethra is the tube that empties the urine out of the body.

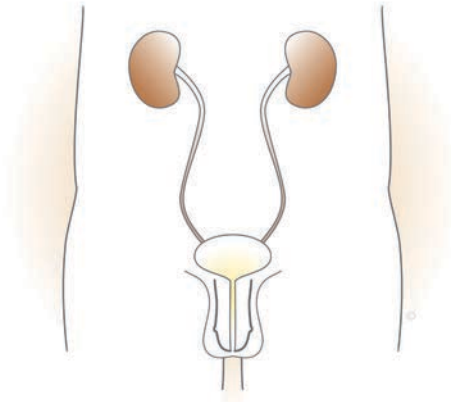

When the bladder is full, the brain sends a signal down the spinal cord to the bladder, causing it to empty. For people with spinal cord issues, the signal from the bladder to the brain gets interrupted, making them unable to empty their bladder.

When people are unable to empty their bladder on their own, they are at risk for urinary tract infections, as well as incontinence or involuntary loss of urine. When urine stays in the bladder and is not emptied, bacteria can grow, causing infections which can lead to illness. Research has shown that self-intermittent catheterization helps reduce urinary tract infections, control urinary leakage (incontinence) and prevent urinary tract damage.

Introduction to Self- Intermittent Catheterization

Intermittent catheterization is the periodic emptying of the bladder by the insertion of a hollow plastic tube (catheter) into the urethra, past the sphincter muscles and into the bladder. Urine then passes out of the bladder through the catheter.

Intermittent catheterization is used when a person is unable to empty his bladder. Medical conditions that often require intermittent catheterization include spinal cord injuries, spina bifida and multiple sclerosis, to name a few.

Intermittent catheterization must be done at regular intervals each day to keep the bladder healthy.

Your healthcare provider has recommended clean self-intermittent catheterization using a coude tip catheter to help empty your bladder and keep your urinary system healthy. A coude tip has a slight angle that makes insertion easier and allows the catheter to move around obstructions such as an enlarged prostate or scar tissue. The Cure Catheter® with coude tip features a blue control stripe to help facilitate proper use.



How Often Should I Catheterize?

Your healthcare provider will let you know how often you will need to perform clean self-intermittent catheterization and the size of the catheter you will need. Normally, self-intermittent catheterization is performed every four hours starting when you wake up in the morning and continues every four hours until bedtime. Some people catheterize on a more frequent schedule. Most people who catheterize do not need catheterization at night. Your healthcare provider will let you know if you need catheterization at night.

Where Do I Get Catheters?

You will be taught clean self-intermittent catheterization by your healthcare provider who will decide the size and style of catheter that you will need. Your healthcare provider will provide a prescription for the catheter supplies. The prescription will be sent to a medical supply company that will provide you with the supplies.

Instructions on Self-Intermittent Catheterization with Coudé Tip Catheter for Men

Self-intermittent catheterization can be performed in a sitting, standing or lying position. Do what is most comfortable for you. If possible, stand in front of the toilet, sit nearby on a chair or wheelchair, or sit on the toilet with your legs spread apart.

1. Before you begin, gather the following materials:

- Cure Catheter® with coudé tip in sealed wrapper
- Water soluble lubricant
- Something to cleanse the area: soap, water and a washcloth, or unscented disposable wipes
- A container to collect the urine if not cathing on or near a toilet
- Hand towel




Male, 16", coudé tip catheter with blue control stripe, polished eyelets, and funnel end shown above.

2. Inspect catheter before use. If catheter or package is damaged do not use.

3. Wash hands thoroughly with soap and water, or if not available, use an antibacterial hand cleaner.

Self-Intermittent Catheterization Instructions for Men *continued*

4. Lay out the equipment so it is within easy reach. Open the water soluble lubricant and catheter package. The Cure Catheter® package opens from both ends for ease of use. Apply a liberal amount of lubricant to the coudé tip and first several inches of the catheter.
5. Position yourself comfortably. Arrange clothing out of the way. If sitting on the toilet, spread your legs apart. If lying down (on towel or underpad) knees can be straight or bent.
6. Hold the penis upright on the sides so that the urethra is not pinched.
7. Retract the foreskin if the penis is uncircumcised. Wash the glans penis from the urethral opening (tip) to the base of the glans with soap, water and washcloth or unscented disposable wipe. Wash with a circular motion from the tip outward. Do not wash back and forth over the urethral opening. Repeat the washing two more times with different areas of the wash cloth or with a new unscented wipe each time.
8. While holding the penis upright with one hand at a 45 to 60-degree angle from the body, slowly insert the lubricated coudé catheter into the penis using your dominant hand with the blue stripe facing upward. Gently slide the catheter, until you meet resistance at the urethral sphincters. Do not pull the catheter in and out when you meet resistance; just apply gentle but firm continuous pressure until the catheter advances into the bladder.
9. Continue to insert the catheter with the blue stripe facing upwards until urine begins to flow. If you are standing in front of the toilet or sitting on it with legs spread apart, the urine can drain into the toilet. If you are lying down, place the funnel end of the catheter into a container to collect the urine. Hold it in place until urine stops flowing.

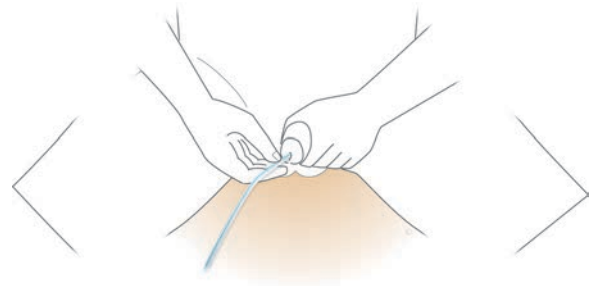

Self-Intermittent Catheterization Instructions for Men continued

10. When urine stops flowing, slowly rotate the catheter between your fingers while withdrawing the catheter. If urine starts to flow again when withdrawing the catheter, stop each time it flows and let the urine drain. Continue to slowly withdraw the catheter until the bladder is empty and the catheter is removed. If you are uncircumcised, make sure you bring the foreskin back over the glans penis when you are finished.

The Cure Catheter® is available in Straight and Coude tip designs with smooth, polished eyelets for maximum user comfort.

Cure Medical catheters are not made with DEHP, DINP, BPA, or natural rubber latex.





Your Clean Self-Intermittent Catheterization Specifications:

- Cure Catheter® with coudé tip size:
 10 FR 12 FR 14 FR 16 FR 18 FR
- Catheterize _____ times a day or every _____ hours.
- Catheterization required at night: Yes No

Healthcare provider contact information:

Name

Telephone Number

Healthcare product supplier contact information:

Name

Telephone Number